

一、 User

01 帳號

i 註冊

METHOD	URL
POST	/api/register

Input

NAME	TYPE	DESCRIPTION
email	String	信箱
password	String	密碼

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗
message	String	成功; 失敗

Example

METHOD		URL	
POST		/api/register	
NAME	TYPE	Value	

email	String	test@gmail.com
password	String	1234

ii 登入

METHOD	URL
POST	/api/auth

Input

NAME	TYPE	DESCRIPTION
email	String	信箱
password	String	密碼

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗; 2=信箱未驗證
token	String (成功才有)	Token
message	String	成功; 失敗

Example

METHOD		URL
POST		/api/auth
NAME	TYPE	Value
email	String	test@gmail.com
password	String	1234

iii 發送驗證碼

METHOD		URL
POST		/api/verification/send

Input

NAME	TYPE	DESCRIPTION
email	String	信箱

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD	URL
--------	-----

POST	/api/verification/send	
NAME	TYPE	Value
email	String	test@gmail.com

iv 檢查驗證碼

METHOD	URL
POST	/api/verification/check

Input

NAME	TYPE	DESCRIPTION
email	String	信箱
code	String	驗證碼

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD	URL
POST	/api/verification/check

NAME	TYPE	Value
email	String	test@gmail.com
code	String	77777

v 忘記密碼

METHOD	URL
POST	/api/password/forgot

Input

NAME	TYPE	DESCRIPTION
email	String	信箱

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD	URL	
POST	/api/password/forgot	
NAME	TYPE	Value

email	String	test@gmail.com
--------------	--------	-----------------------

vi 重設密碼

METHOD		URL
POST	/api/password/reset	
HEADER	VALUE	DESCRIPTION
Authorize	sessionid=\$token	Token

Input

NAME	TYPE	DESCRIPTION
password	String	密碼

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD		URL
POST	/api/password/reset	
HEADER	VALUE	DESCRIPTION

Authorization	Bearer \$token	Token
NAME	TYPE	Value
password	String	密碼

vii 註冊確認

METHOD	URL
GET	/api/register/check

Input

NAME	TYPE	DESCRIPTION
email	String	信箱

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD		URL	
GET	/api/register/check		
NAME	TYPE	Value	

email	String	test@gmail.com
-------	--------	----------------

02 其他

i 最新消息

METHOD	URL	
GET	/api/news	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token

Input

NAME	TYPE	DESCRIPTION
------	------	-------------

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗
news	Object	最新消息
<pre>[{ "id": 1, // int</pre>		


```
"member_id": 1, // int

"group": 1, // int

"title": "test", // str

"message": "test", // str

"pushed_at": "2023-8-29 16:30:30", // str

"created_at": "2023-8-29 16:30:30", // str

"updated_at": "2023-8-29 16:30:30" // str

},

...

]
```

ii 分享

METHOD		URL	
POST		/api/share	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
type	Int	0: 血壓; 1: 體重; 2: 血糖; 3: 飲食;
id	Int	紀錄 ID
relation_type	Int	1: 親友; 2: 糖友

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD		URL	
POST		/api/share	

HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token
NAME	TYPE	Value
type	Int	0
id	Int	1
relation_type	Int	1

iii 查看分享 (找不此 API 觸發方式)

METHOD		URL	
GET	/api/share/<type>		
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	
說明			
將根據指定的關係，回傳所有記錄的資料後，把所有資料個別整理成血壓、體重、血糖、減肥日記，用 JSON 格式回傳。			

Input

NAME	TYPE	DESCRIPTION
type	String	0: 醫師; 1: 親友; 2: 糖友

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗
records	Object	紀錄

二、 Body

01 個人資訊

i 個人資訊設定

METHOD		URL	
PATCH		/api/user	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	
說明			
設定 UserProfile 與 Userset 資料			

Input

NAME	TYPE	DESCRIPTION
name	String	暱稱
birthday	String	生日
height	Double	身高
weight	String	體重
phone	String	電話
email	String	信箱

gender	Bool	0: 男; 1: 女
fcm_id	String	Fcm id
address	String	地址

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD	URL	
PATCH	/api/user	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token
NAME	TYPE	Value
name	String	王小明
birthday	String	2023-08-29
height	Double	187
weight	String	87
phone	String	0987987987
email	String	test@gmail.com

gender	Bool	true
fcm_id	String	1234
address	String	台中市北區三民路三段 129 號

ii 個人資訊

METHOD		URL	
GET	/api/user		
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	
說明			
讀取 UserProfile、Userset、Default、Setting 的資料			

Input

NAME	TYPE	DESCRIPTION
------	------	-------------

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗
user	Object	使用者資訊

```
{  
  
  "id": 1, // int  
  
  "name": "王小明", // str  
  
  "account": "test@gmail.ocm", // str  
  
  "email": "test@gmail.ocm", // str  
  
  "phone": "0987987987", // str  
  
  "fb_id": "未設置", // str  
  
  "status": "VIP", // str  
  
  "group": "0", // str  
  
  "birthday": "2023-08-29", // str  
  
  "height": , 187, // double  
  
  "weight": "87", // str  
  
  "gender": 0, // int  
  
  "address": "台中市北區三民路三段 129 號", // str  
  
  "unread_records": [ // list  
  
    0, // int  
  
    0, // int  
  
    0 // int  
  
  ],
```


"verified": 0, // int

"privacy_policy": 1, // int

"must_change_password": 1, // int

"fcm_id": "", // str

"login_times": 0, // int

"created_at": "2023-08-23 16:51:14" // str

"updated_at": "2023-08-23 16:51:14" // str

"default": {

"id": 1, // int

"user_id": 1, // int

"sugar_delta_max": null, // double

"sugar_delta_min": null, // double

"sugar_morning_max": null, // double

"sugar_morning_min": null, // double

"sugar_evening_max": null, // double

"sugar_evening_min": null, // double

"sugar_before_max": null, // double

"sugar_before_min": null, // double

"sugar_after_max": null, // double

```
"sugar_after_min": null, // double

"systolic_max": null, // int

"systolic_min": null, // int

"diastolic_max": null, // int

"diastolic_min": null, // int

"pulse_max": null, // int

"pulse_min": null, // int

"weight_max": null, // double

"weight_min": null, // double

"bmi_max": null, // double

"bmi_min": null, // double

"body_fat_max": null, // double

"body_fat_min": null, // double

"created_at": "2023-08-23 16:51:14", // str

"updated_at": "2023-08-23 16:51:14" // str

},

"setting": {

  "id": 1, // int

  "user_id": 1, // int
```

```
"after_recording": 0, // int

"no_recording_for_a_day": 0, // int

"over_max_or_under_min": 0, // int

"after_meal": 0, // int

"unit_of_sugar": 0, // int

"unit_of_weight": 0, // int

"unit_of_height": 0, // int

"created_at": "2023-02-03 08:17:17", // str

"updated_at": "2023-02-03 08:17:17" // str

},

"vip": {

  "id": 1, // int

  "user_id": 1, // int

  "level": 0, // int

  "remark": 0.0, // double

  "started_at": "2023-02-03 08:17:17", // str

  "ended_at": "2023-02-03 08:17:17", // str

  "created_at": "2023-02-03 08:17:17", // str

  "updated_at": "2023-02-03 08:17:17", // str
```

```
}  
  
}
```

iii 個人預設值

METHOD		URL	
PATCH		/api/user/default	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	
說明			
將預設值更新到 Default 資料表			

Input

NAME	TYPE	DESCRIPTION
sugar_delta_max	Double	預設值
sugar_delta_min	Double	預設值
sugar_morning_max	Double	預設值
sugar_morning_min	Double	預設值
sugar_evening_max	Double	預設值
sugar_evening_min	Double	預設值
sugar_before_max	Double	預設值

sugar_before_min	Double	預設值
sugar_after_max	Double	預設值
sugar_after_min	Double	預設值
systolic_max	Int	預設值
systolic_min	Int	預設值
diastolic_max	Int	預設值
diastolic_min	Int	預設值
pulse_max	Int	預設值
pulse_min	Int	預設值
weight_max	Double	預設值
weight_min	Double	預設值
bmi_max	Double	預設值
bmi_min	Double	預設值
body_fat_max	Double	預設值
body_fat_min	Double	預設值

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD		URL	
PATCH	/api/user/default		
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	
NAME	TYPE	Value	
sugar_delta_max	Double	10.0	
sugar_delta_min	Double	1.0	
sugar_morning_max	Double	10.0	
sugar_morning_min	Double	1.0	
sugar_evening_max	Double	10.0	
sugar_evening_min	Double	1.0	
sugar_before_max	Double	10.0	
sugar_before_min	Double	1.0	
sugar_after_max	Double	10.0	
sugar_after_min	Double	1.0	
systolic_max	Int	10	
systolic_min	Int	1	
diastolic_max	Int	10	

diastolic_min	Int	1
pulse_max	Int	10
pulse_min	Int	1
weight_max	Double	10.0
weight_min	Double	1.0
bmi_max	Double	10.0
bmi_min	Double	1.0
body_fat_max	Double	10.0
body_fat_min	Double	1.0

iv 個人設定

METHOD		URL	
PATCH		/api/user/setting	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
after_recording	Bool	設定値

no_recording_for_a_day	Bool	設定値
over_max_or_under_min	Bool	設定値
after_meal	Bool	設定値
unit_of_sugar	Bool	設定値
unit_of_weight	Bool	設定値
unit_of_height	Bool	設定値

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD	URL	
PATCH	/api/user/setting	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token
NAME	TYPE	Value
after_recording	Bool	false

no_recording_for_a_day	Bool	false
over_max_or_under_min	Bool	true
after_meal	Bool	true
unit_of_sugar	Bool	true
unit_of_weight	Bool	true
unit_of_height	Bool	true

02 測量上傳

i 上傳血壓測量結果

METHOD		URL
POST	/api/user/blood/pressure	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token

Input

NAME	TYPE	DESCRIPTION
Systolic	Int	收縮壓

Diastolic	Int	舒張壓
Pulse	Int	心跳
Recorded_at	String	紀錄時間

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗
records	String	成功; 失敗

Example

METHOD	URL	
POST	/api/user/blood/pressure	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token
NAME	TYPE	Value
Systolic	Int	100
Diastolic	Int	100
Pulse	Int	60
Recorded_at	String	2023-02-03 08:17:17

ii 上傳體重測量結果

METHOD		URL	
POST		/api/user/weight	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
weight	Double	體重
body_fat	Double	體脂
bmi	Double	BMI
recorded_at	String	紀錄時間

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD	URL
POST	/api/user/weight

HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token
NAME	TYPE	Value
weight	Double	87.0
body_fat	Double	23.0
bmi	Double	23.0
recorded_at	String	2023-02-03 08:17:17

iii 上傳血糖

METHOD	URL	
POST	/api/user/blood/sugar	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token

Input

NAME	TYPE	DESCRIPTION
sugar	Double	血糖
timeperiod	Int	時段
recorded_at	String	紀錄時間

drug	Int	有使用藥物
exercise	Int	有運動

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD	URL	
POST	/api/user/blood/sugar	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token
NAME	TYPE	Value
sugar	Double	123.0
timeperiod	Int	1
recorded_at	String	2023-02-03 08:17:17
drug	Int	1
exercise	Int	1

iv 上一筆紀錄資訊

METHOD		URL	
POST		/api/user/records	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
diet	Int	時段

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗
blood_sugars	Object	
{ "sugar": 0.0 // double }		
blood_pressures	Object	
{		

<pre> "systolic": 345, // int "diastolic": 345, // int "pulse": 345 // int } </pre>		
weights	Object	
<pre> { "weight": 0.0 // double } </pre>		

Example

METHOD		URL	
POST		/api/user/records	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	
NAME	TYPE	Value	
diet	Int	時段	

03 日記

i 日記列表資料

METHOD		URL	
GET		/api/user/diary	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
date	String	日期

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗
diary	Object	日記

```
[  
  
  {  
  
    "id": 1, // int  
  
    "user_id": 1, // int
```



```
"systolic": 140, // int

"diastolic": 139, // int

"pulse": 117, // int

"weight": 0.0, // double

"body_fat": 0.0, // double

"bmi": 0.0, // double

"sugar": 0.0, // double

"exercise": 0, // int

"drug": 0, // int

"timeperiod": 0, // int

"description": "", // str

"meal": 3, // int

"tag": [

    "name": ["abc"], // list[str]

    "message": "" // str

],

"image": [""], // list[str]

"location": {

    "lat": "", // str
```

```
    "lng": "" // str

  },

  "reply": "", // str

  "recorded_at": "2023-09-21 00:06:14", // str

  "type": "blood_pressure" // str

}

]
```

Example

METHOD		URL	
GET		/api/user/diary/	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	
NAME	TYPE	Value	
date	String	2023-09-21	

ii 飲食日記

METHOD		URL	
POST		/api/user/diet	

HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token

Input

NAME	TYPE	DESCRIPTION
description	String	描述
meal	Int	時段
tag[]	List [String]	Tags
image	Int	照片數量
lat	Double	緯度
lng	Double	經度
recorded_at	String	使用者的記錄時間

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗
image_url	String	圖片網址

Example

METHOD	URL
POST	/api/user/diet

HEADER	VALUE	HEADER
Authorization	Bearer \$token	Authorization
NAME	TYPE	Value
description	String	good
meal	Int	0
tag[]	List [String]	dinner
image	Int	1
lat	Double	120.5
lng	Double	323.5
recorded_at	String	2023-09-21 00:06:14

iii 刪除日記記錄

METHOD	URL	
DELETE	/api/user/records	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token

Input

NAME	TYPE	DESCRIPTION
------	------	-------------

deleteObject	List [String] (Optional)	IDs
---------------------	---------------------------------	-----

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD	URL	
DELETE	/api/user/records	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token
NAME	TYPE	Value
blood_sugars[]	Param	1

04 糖化血色素

i 糖化血色素

METHOD	URL
GET	/api/user/a1c

HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token

Input

NAME	TYPE	DESCRIPTION
------	------	-------------

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗
a1cs	List [Object]	糖化血色素

[
{
"id": 1, // int
"user_id": "", // int
"a1c": "10", // str
"recorded_at": "2023-02-07 04:54:43", // str
"created_at": "2023-02-07 04:54:44", // str
"updated_at": "2023-02-07 04:54:44" // str
},
...

1

ii 送糖化血色素

METHOD		URL	
POST		/api/user/a1c	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
A1c	Int	糖化血色素
Recorded_at	String	紀錄時間

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD	URL
POST	/api/user/a1c

HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token
NAME	TYPE	Value
A1c	Int	10
Recorded_at	String	2023-02-07 04:54:43

iii 刪除糖化血色素

METHOD		URL	
DELETE		/api/user/a1c	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
ids[]	Dict [List <Int>]	IDs

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD		URL	
DELETE		/api/user/a1c	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	
NAME	TYPE	Value	

ids[]	Dict [List <Int>]	{ "ids[]": 52] }
--------------	--------------------	---

05 就醫資訊

i 就醫資訊

METHOD		URL	
GET		/api/user/medical	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
------	------	-------------

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

medical_info	Object	就醫資訊
<pre> { "id": 1, // int "user_id": 1, // int "diabetes_type": 1, // int "oad": 0, // int "insulin": 0, // int "anti_hypertensives": 0, // int "created_at": "", // str "updated_at": "" // str }</pre>		

ii 更新就醫資訊

METHOD		URL	
PATCH		/api/user/medical	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
diabetes_type	Int	0~4 無/糖尿病前期/第一型/ 第二型/妊娠
oad	Bool	糖尿病口服藥
insulin	Bool	胰島素
anti_hypertensives	Bool	高血壓藥

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD	URL	
PATCH	/api/user/medical	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token
NAME	TYPE	Value
diabetes_type	Int	0
oad	Bool	true

insulin	Bool	true
anti_hypertensives	Bool	true

06 藥物資訊

i 藥物資訊

METHOD		URL
GET		/api/user/drug-used
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token

Input

NAME	TYPE	DESCRIPTION
------	------	-------------

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗
drug_useds	Object	藥物資訊
<pre>[{</pre>		

```
"id": 0, // int

"user_id": 0, // int

"type": 1, // int

"name": "", // str

"recorded_at": "", // str

"updated_at": "", // str

"created_at": "", // str

},

...

]
```

ii 上傳藥物資訊

METHOD		URL	
POST		/api/user/drug-used	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
type	Int	0: 糖尿病藥物

		1: 高血壓藥物
name	String	藥物的名稱
recorded_at	String	紀錄時間

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD	URL	
POST	/api/user/drug-used	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token
NAME	TYPE	Value
type	Int	0
name	String	test
recorded_at	String	2023-02-07 04:54:43

iii 刪除藥物資訊

METHOD		URL	
DELETE		/api/user/drug-used	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
ids[]	Dict [String, List<Int>]	IDs

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD		URL	
DELETE		/api/user/drug-used	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

NAME	TYPE	Value
ids[]	Dict [String, List<Int>]	{ "ids[]":[52] }

07 關懷資訊

i 獲取關懷諮詢

METHOD	URL	
GET	/api/user/care	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token

Input

NAME	TYPE	DESCRIPTION
------	------	-------------

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗
cares	Object	關懷資訊

```
[
  {
    "id": 0, // int
    "user_id": 1, // int
    "member_id": 1, // int
    "reply_id": null, // int reply_id=null 為罐頭訊息
    "message": "", // str
    "created_at": "2017-11-30 15:46:11", // str
    "updated_at": "2017-11-30 15:46:11" // str
  },
  ...
]
```

ii 發送關懷諮詢

METHOD		URL	
POST		/api/user/care	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
Message	String	關懷訊息

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD	URL	
POST	/api/user/care	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token
NAME	TYPE	Value
Message	String	hi

08 其他

i 更新 Badge

METHOD	URL
PUT	/api/user/badge

HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token

Input

NAME	TYPE	DESCRIPTION
Badge	Int	Badge

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD	URL	
PUT	/api/user/badge	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token
NAME	TYPE	Value
Badge	Int	0

三、 Friend

01 控糖團好友

i 獲取控糖團邀請碼

METHOD		URL	
GET		/api/friend/code	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
------	------	-------------

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗
invite_code	String	控糖團邀請碼

ii 控糖團列表

METHOD	URL
--------	-----

GET	/api/friend/list	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token

Input

NAME	TYPE	DESCRIPTION
------	------	-------------

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗
friends	Object	好友列表

```
[
  {
    "id": 1,  // int
    "name": "",  // str
    "relation_type": 0  // int
  }
]
```

iii 獲取控糖團邀請

METHOD		URL	
GET		/api/friend/requests	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
------	------	-------------

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗
requests	Object	控糖團邀請

```
[
  {
    "id": 1, // int
    "user_id": 2, // int
    "relation_id": 1, // int
    "type": 0, // int
  }
]
```

```
"status": 0, // int

"read": 0, // int

"created_at": "2017-10-20 15:46:16", // str

"updated_at": "2017-10-20 15:49:20", // str

"user": {

    "id": 2, // int

    "name": null, // str

    "account": "fb_2" // str

}

}

]
```

iv 送出控糖團邀請

METHOD		URL	
POST	/api/friend/send		
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
type	Int	0: 醫師團; 1: 親友團; 2: 糖友團
invite_code	String	邀請碼

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=邀請碼無效, 2=已經成為好友
message	String	成功; 失敗

Example

METHOD	URL	
POST	/api/friend/send	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token
NAME	TYPE	Value
type	Int	0
invite_code	String	55688

v 接受控糖團邀請

METHOD		URL	
GET		/api/friend/<inviteid>/accept	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
inviteid	Int	ID

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD		URL	
GET		/api/friend/<inviteid>/accept	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	
NAME	TYPE	Value	

inviteld	Int	55688
-----------------	-----	-------

vi 拒絕控糖團邀請

METHOD		URL	
GET		/api/friend/<inviteId>/refuse	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	

Input

NAME	TYPE	DESCRIPTION
inviteId	Int	ID

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗

Example

METHOD		URL	
GET		/api/friend/<inviteId>/refuse	
HEADER	VALUE	DESCRIPTION	
Authorization	Bearer \$token	Token	
NAME	TYPE	Value	

inviteId	Int	55688
-----------------	-----	-------

vii 删除更多好友 OK

METHOD		URL
DELETE	/api/friend/remove	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token

Input

NAME	TYPE	DESCRIPTION
ids[]	Dict[Str, Int]	ID

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失败;
message	String	成功; 失败

Example

METHOD		URL
DELETE	/api/friend/remove	
HEADER	VALUE	DESCRIPTION

Authorization	Bearer \$token	Token
NAME	TYPE	Value
ids[]	Dict[Str, Int]	<code>{"ids[]": [52]}</code>

viii控糖團結果 OK

METHOD	URL	
GET	/api/friend/results	
HEADER	VALUE	DESCRIPTION
Authorization	Bearer \$token	Token

Input

NAME	TYPE	DESCRIPTION
------	------	-------------

Output

NAME	TYPE	DESCRIPTION
status	String	0=成功; 1=失敗;
message	String	成功; 失敗
results	List [Object]	
<pre>[{</pre>		

```
"id": 1, // int

"user_id": 3, // int

"relation_id": 4, // int

"type": 1, // int 0: 醫師 1: 親友 2: 糖友

"status": 1, // int 1: 接受 2: 拒絕

"read": 1, // int

"created_at": "2017-10-30 19:08:26", // str

"updated_at": "2017-10-30 19:18:18", // str

"relation": {

    "id": 4, // int

    "name": null, // str

    "account": "fb_1" // str


}

},


...

]
```


四、 流程圖(參考就好)

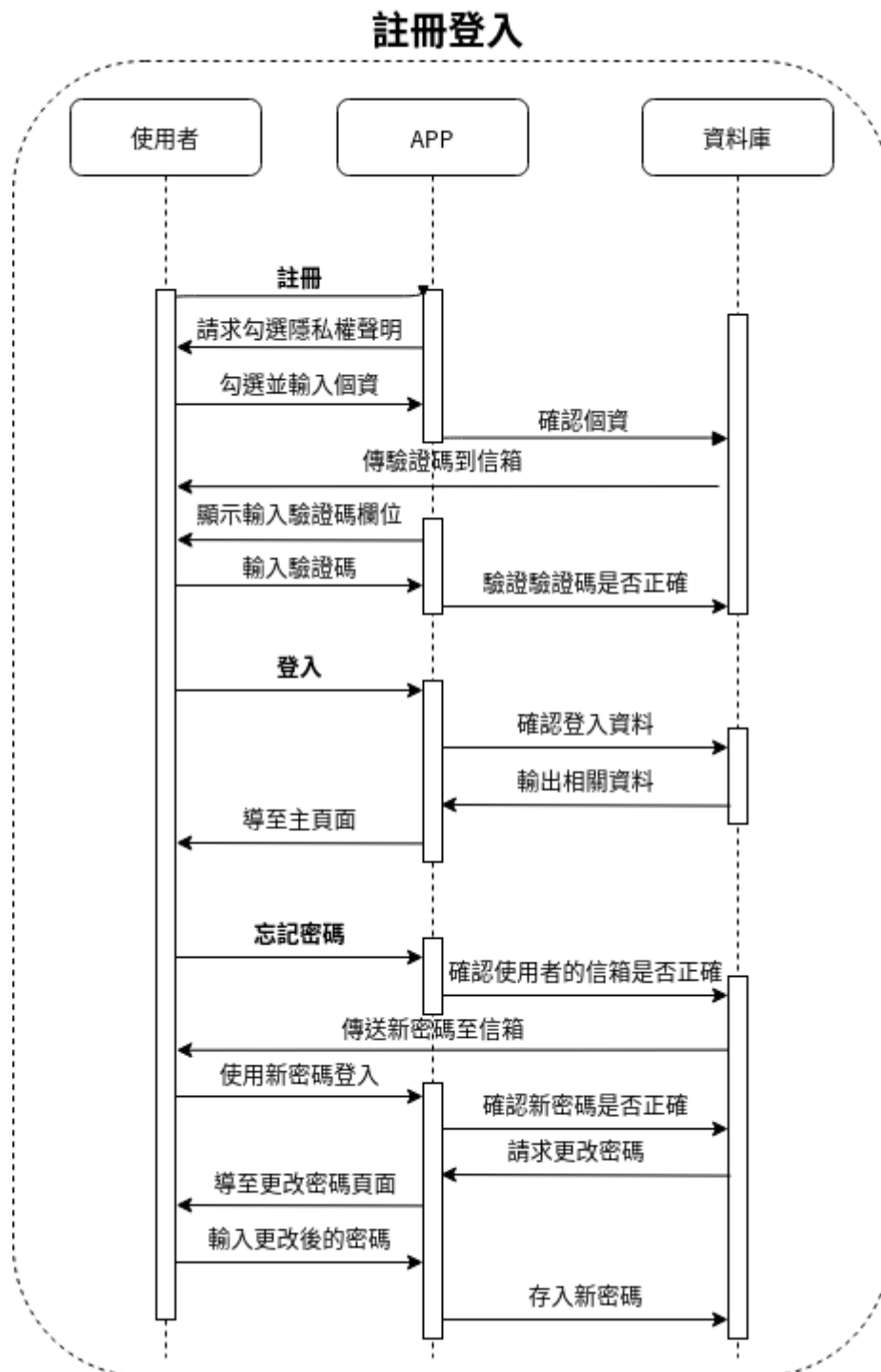


個人資訊

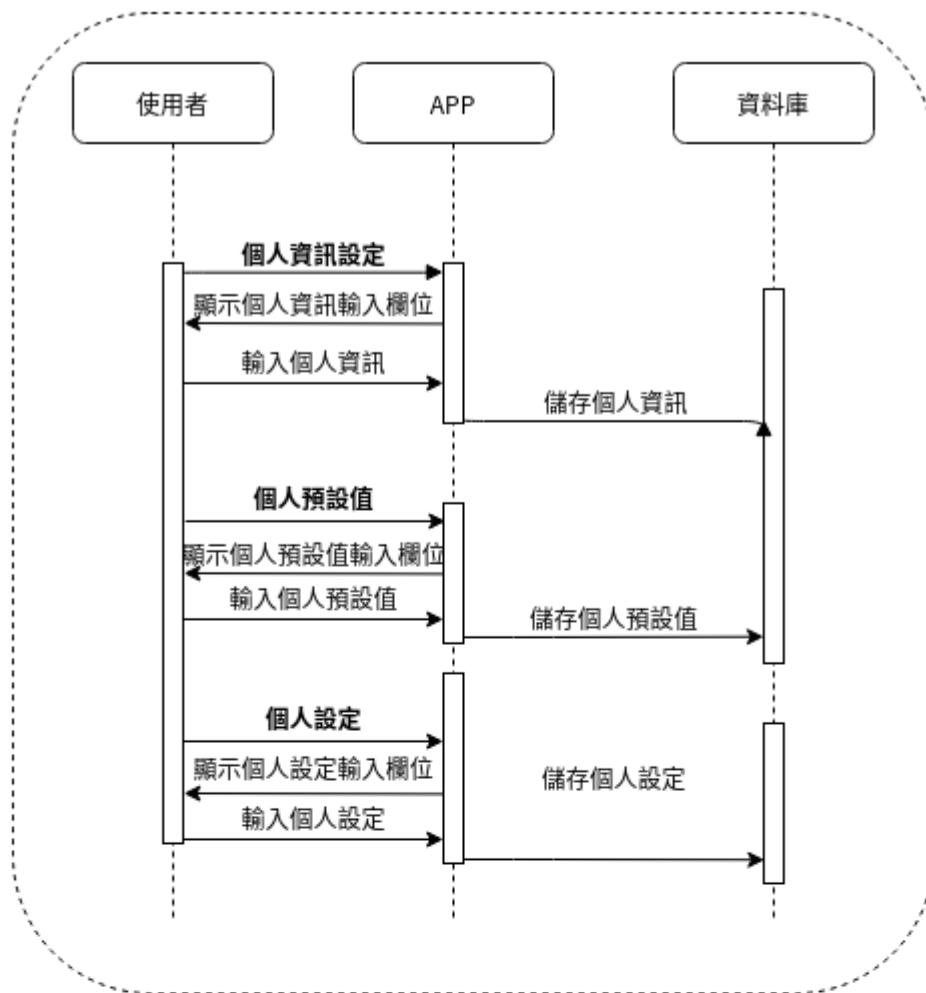


糖化血色素

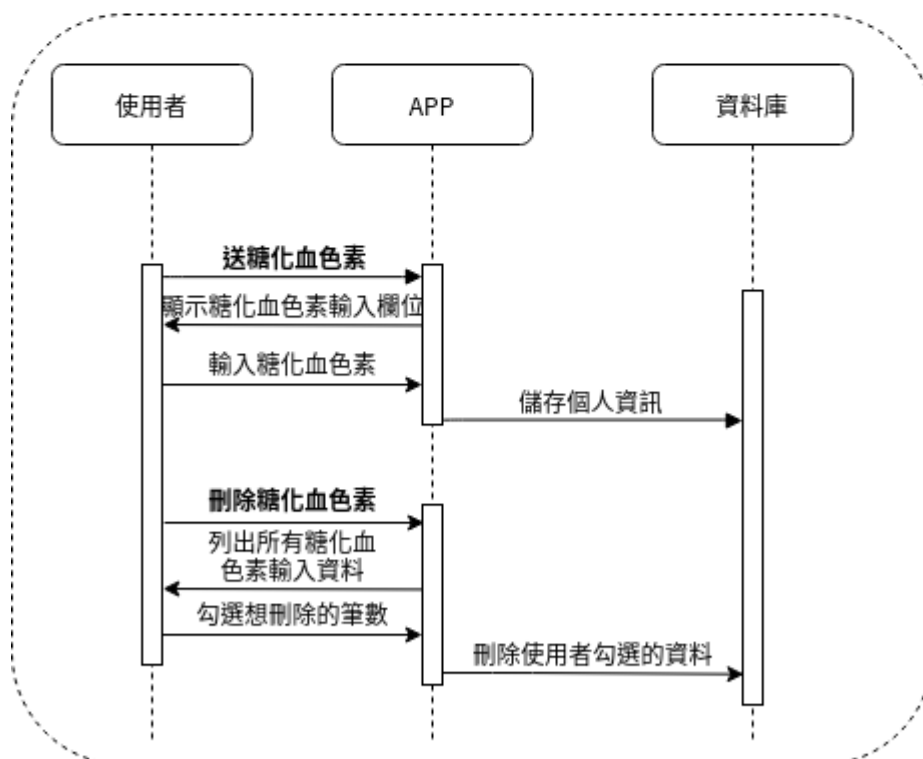


就醫、藥物資訊

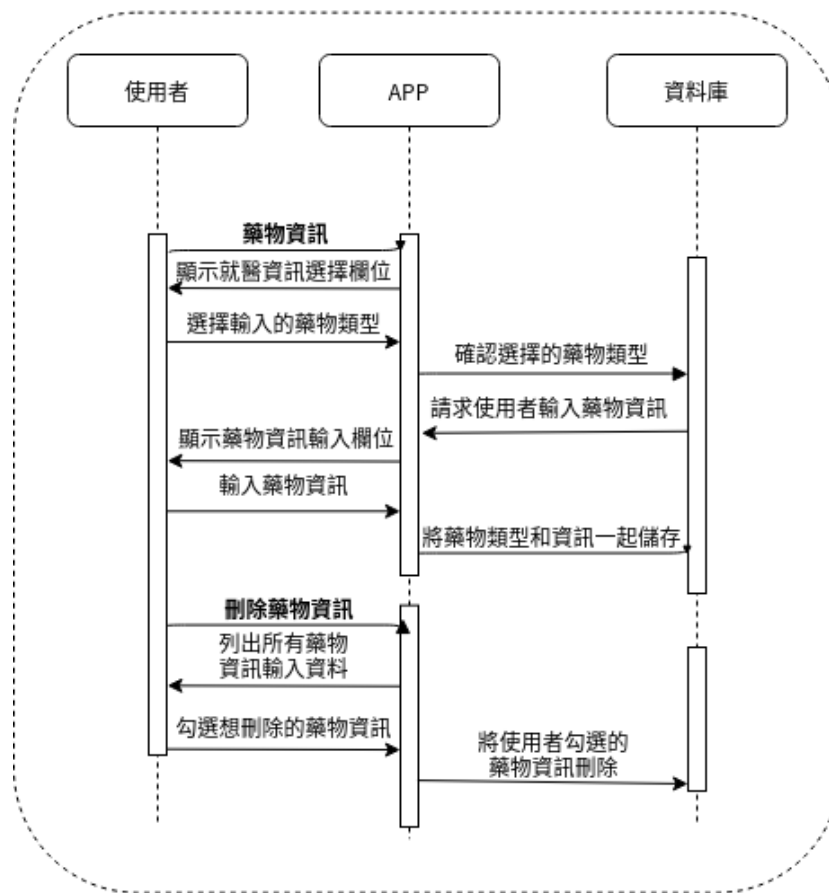

日記

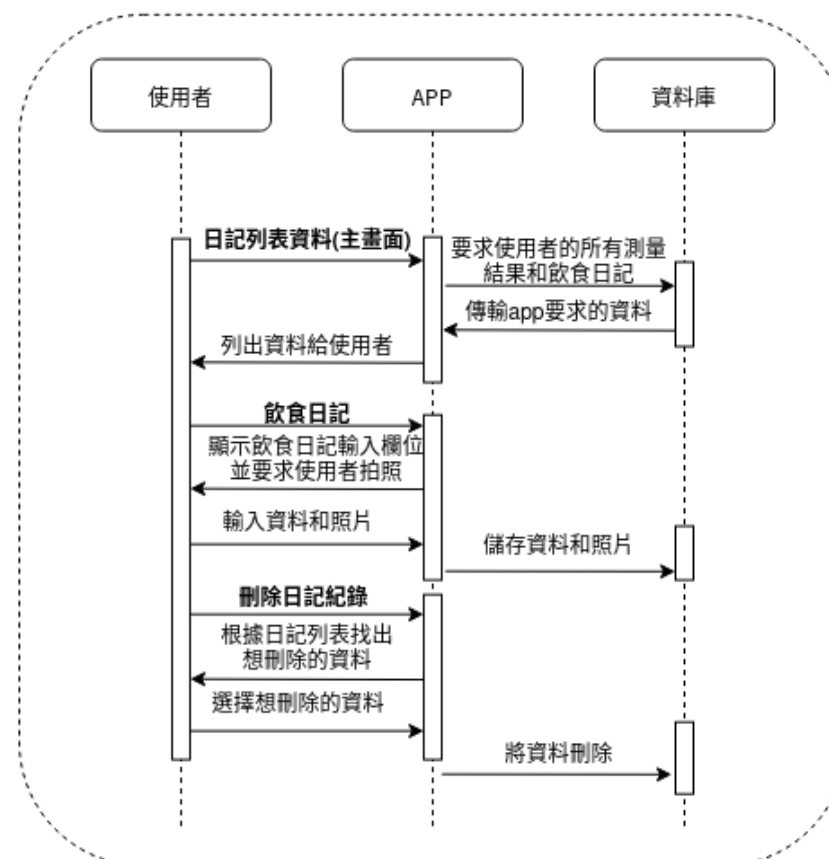

控糖團

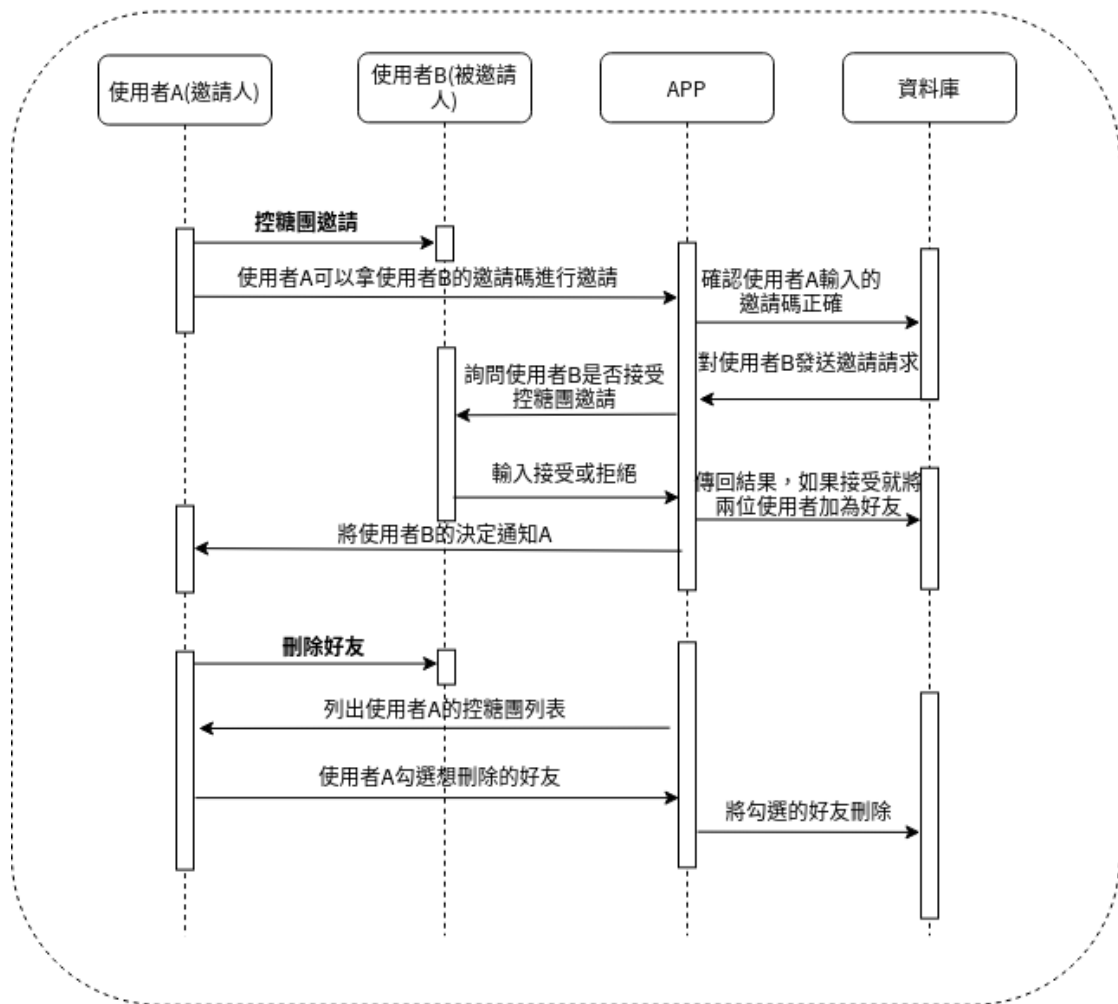

測量上傳

